

BLICKMAN



IV Stands

Optimal continuous fluid delivery

Includes:

- IV Stands
- IV Stand Accessories
- IV Irrigation Safety Systems




THE BLICKMAN DIFFERENCE

IV Stands

- Our IV poles are optimal for maintaining continuous fluid delivery for every patient, available in a wide range of adjustable heights with a variety of easy locking options that can be operated with one hand.
- Our IV poles are all available in high quality chrome and stainless steel, and set on 4, 5, or 6 leg bases with casters, no projections, and great balance/tip resistance. Some bases are gravity weighted to help compensate for heavier loads. All bases have the ability to spin 360° in tight spaces.
- Most of our models have 2-, 4-, 6- and 8-(rake) hook adapters available. All come with colored secure grip tips.

Key Product Features

*See page 9 for more details on all optional features

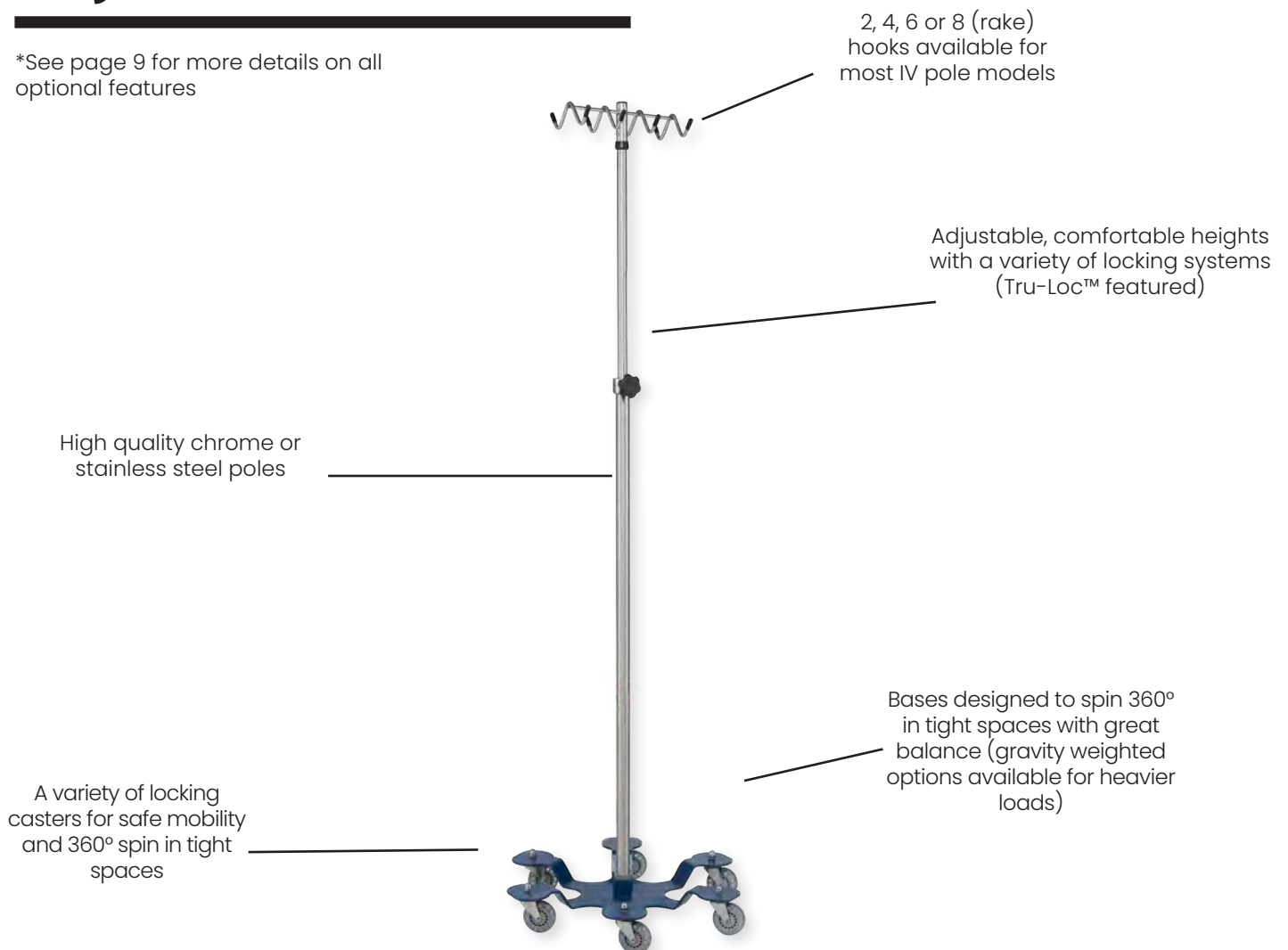


IV Stands

4 Leg

4 LEG CHROME IV STANDS

Description: We offer two 4-leg chrome IV stand models, one with a standard base, and one with a heavyweight base. Both base types include Wall Saver™ tips to protect against wall scuff marks. They both move on 2" plastic twin wheel hooded swivel casters and they both come standard with two hooks. A four hook removable adapter is available for model 1305.

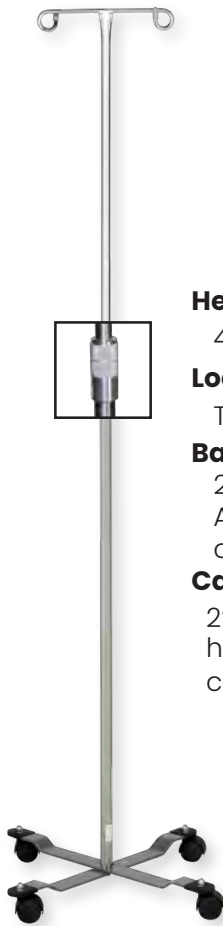

Height adjustment:

47 1/2" - 85 1/2"

Lock:

Twist

Base:

22" diameter
Apoxy coated
chrome

Casters

2" plastic twin wheel
hooded swivel
casters



Height adjustment

51 1/2" - 93"

Lock:

Tru-loc™

Base:

22" diameter
Chrome plated
heavy bar steel

Casters

2" plastic twin wheel
hooded swivel
casters


Name	Standard Base	Heavyweight Base
# of Hooks	2 Welded Hooks	2 Removable Hooks
Model #	1305	1310
Item #	0561305000	0561310000
Shipping Weight	13 lbs.	14 lbs.

IV Stands

4 Leg

4 LEG STAINLESS STEEL IV STANDS

Description: We offer two types of 4-leg stainless steel IV stands, one with a standard base, and one with a heavyweight base. They come standard with either 2 or 4 removable hooks. Contact customer service for 6- or 8-(rake) hooks.

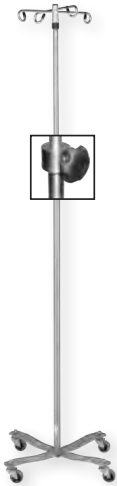

Name	Standard Base
# of Hooks	2 Hooks
Model #	1410SS
Item #	0571410000
Shipping Weight	13 lbs.
# of Hooks	4 Hooks
Model #	1410SS-4
Item #	0571410400
Shipping Weight	14 lbs.

Height Adjustment: 51 ³/₄" - 94"
 Mechanism: Tru-Loc™
 Base Material: Stainless steel with Wall Savers™, 22" diameter
 Casters: 2" hard rubber swivel casters



Name	Heavyweight Base	
Lock Type	Twist Lock	Thumb Control
# of Hooks	2 Hooks	
Model #	8890SS	8889SS
Item #	0528890000	0518889000
Shipping Weight	24 lbs.*	24 lbs.*
# of Hooks	4 Hooks	
Model #	8890SS-4	8889SS-4
Item #	0528890400	0518889400
Shipping Weight	24 lbs.*	24 lbs.*

Height Adjustment: 66" - 100"
 Mechanism: Twist Lock, Thumb Control, Foot Operated
 Base Material: Heavyweight stainless steel flat bar, 22" diameter
 Casters: 2" roller blade swivel casters

IV Stands

5 Leg


5 LEG IV STANDS

Description: Our 5-leg IV stands come in either chrome (standard & heavy duty) or stainless steel poles and stand on a 24" diameter, polished die cast aluminum base with Wall Saver™ caps. Contact customer service for 6- or 8-(rake) hooks on certain models.



Name	Chrome Pole	
	Pole	
	Standard (1")	Heavy Duty (1 1/4")
# of Hooks	2 Hooks	4 Hooks
Caster Heights	2 1/2" Plastic	3" Rubber
Height Adjustments	52 1/2" - 93 1/2"	56" - 100"
Model #	1315-2	1350-4
Item #	0561315000	0561350400
Shipping Weight	16 lbs.*	17 lbs.*
	4 Hooks	6 Hooks
Caster Heights	2 1/2" Plastic	3" Rubber
Height Adjustments	52 1/2" - 93 1/2"	56" - 100"
Model #	1315-4	1350-6
Item #	0561315400	0561350600
Shipping Weight	16 lbs.*	17 lbs.*


Height Adjustment Mechanism: Tru-Loc™
 Base Material: 5 Leg, Polished die cast aluminum, 24" diameter
 Casters: 3" hard rubber OR 2 1/2" twin wheel swivel casters

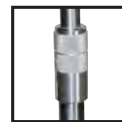

Name

Stainless Steel Pole

Lock Type	Tru-Loc™	Twist Lock	Thumb Control
# of Hooks		2 Hooks	
Caster Heights	2" Rubber	3" Rubber	3" Rubber
Height Adjustments	52" - 94"	74" - 110"	67" - 98"
Model #	1415SS	7792SS	7794SS
Item #	0531415000	0537792000	0517794000
Shipping Weight	16 lbs.*	15 lbs.*	16 lbs.*
# of Hooks		4 Hooks	
Caster Heights	2" Rubber	3" Rubber	3" Rubber
Height Adjustments	52" - 94"	74" - 110"	67" - 98"
Model #	1415SS-4	7792SS-4	7794SS-4
Item #	0531415400	0537792400	0517794400
Shipping Weight	16 lbs.*	16 lbs.*	17 lbs.*

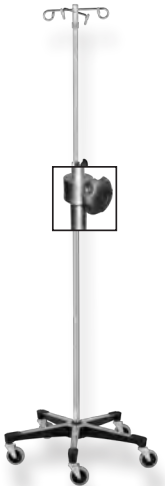

Height Adjustment Mechanism: Tru-Loc™, Twist Lock, Thumb Control
 Base Material: 5 Leg, Polished die cast aluminum, 24" diameter
 Casters: 2" or 3" hard rubber swivel casters


*Shipping weight with an asterisk is the combined shipping weight of separate parts

MR option available, see MR section for more details



IV Stands

5 Leg

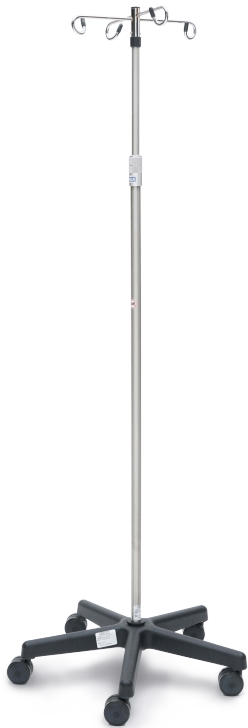

5 LEG WASHABLE IV STAND

Description: Our Hackensack 5-leg IV stand is cart washable as it was built specifically to withstand high temperatures and cleaning fluids with weep holes in the pole and in the base for fluid drainage. This model comes with a stainless steel pole with a twist lock height adjustment mechanism. Choose from 2-, 4-, and 8-(rake) removable hooks.

Features Include:

- ✓ 22 7/8" diameter, 5 leg black composite base
- ✓ 2 1/2" twin wheel hooded swivel casters
- ✓ 1" outside diameter base pole, and 3/4" outside diameter extension pole
- ✓ Height extends from 74" to 110" off the ground



Name

Washable IV Stand

Lock Type

Twist Lock

of Hooks

2 Hooks



Model #

7795SS-2

Item #

0539A02302

Shipping Weight

16 lbs.*

of Hooks

4 Hooks



Model #

7795SS-4

Item #

0539A02304

Shipping Weight

16 lbs.*

of Hooks

8 Hooks



Model #

7795SS-RH

Item #

0539A023RH

Shipping Weight

16 lbs.*

IV Stands

6 Leg

6 LEG IV STANDS



Description: Our 6-leg IV stands are designed for hospital applications requiring heavier loads by using low center of gravity bases. They come in either a chrome or a stainless steel pole with a Tru-Loc™ height adjustment mechanism that allows the pole to reach a minimum of 52" high off the ground and a maximum of 90" off the ground. Our 6 leg IV stands come standard with 4-, 6-, and 8-(rake) removable hooks.

Features Include:

- ✓ 23 3/8" diameter base with low center of gravity for improved stability
- ✓ Heavy duty 3" roller blade swivel casters for easy maneuverability
- ✓ Steel powder coated finished base for durability (many colors available)
- ✓ Accommodates up to 4 infusion pumps and up to 8 IV solution bags while preventing tubing line entanglement



Powder Coated Base

Stainless Steel Base

Name	Chrome Pole	SS Pole	SS Pole
# of Hooks	4 Hooks		
Model #	1370-4	1370SS-4	1371SS-4
Item #	0561370400	0541370400	0541371400
Shipping Weight	29 lbs.*	29 lbs.*	26 lbs.*
# of Hooks	6 Hooks		
Model #	1370-6	1370SS-6	1371SS-6
Item #	0561370600	0541370600	0541371600
Shipping Weight	29 lbs.*	29 lbs.*	26 lbs.*
# of Hooks	8 (rake) Hooks		
Model #	1370-RH	1370SS-RH	1371SS-RH
Item #	0561370RH1	05413700RH	05413710RH
Shipping Weight	29 lbs.*	29 lbs.*	26 lbs.*



Custom Powder Coated Colors Available Upon Request

Please email us at info@blickman.com

IV Stands

Options

IV Hook Configurations



Part #	Description	Shipping Weight
780063700L	2 hooks	2 lbs.
780063800L	4 hooks	2 lbs.
780066900L	6 hooks	2 lbs.
7802808750	8 (rake) hooks	6 lbs.

Lock Options



Thumb Lock

Tru-Lock™

Twist Lock

Base Options

4 Leg Chrome Bases



Standard



Heavyweight

4 Leg Stainless Steel Bases



Standard



Heavyweight

5 Leg Bases



Standard



Washable

6 Leg Bases



Powder Coated



Stainless Steel


IV Stands

Accessories

IV POLE ATTACHMENTS



Wire Basket

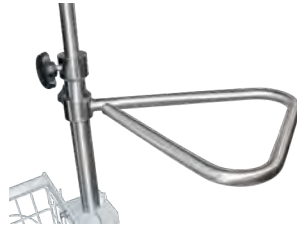

Guidance Handle



Power Strip



Accessory Clamp



Drainage Bag Hook



Pump Support Tray



Multi-pump Support



Tank holder

Part #	Description	Shipping Weight
3501389001	Wire Basket	9 lbs.
3501381000	Guidance Handle (MR Option Available)	4 lbs.
9101107600	Power Strip	4 lbs.
7801380000	Accessory Clamp	4 lbs.
7801382000	Drainage Bag Hook	1 lbs.
7801384000	Pump Support Tray	2 lbs.
7801385000	Multi-pump Support	2 lbs.
7801386000	Tank Holder	5 lbs.

MR options available where indicated, see MR section for more details.



IV Stands

Construction Specifications

Construction Specifications

Universal Dimensions	Height Range: 47" - 110" Base Diameter Range: 22" - 24"
Materials	304 stainless steel Chrome
Gauge	Stainless Steel IV Stands: Base Pole: 1" Diameter, 16GA; Extension Pole: 1" Diameter, 18GA Exception: 1370/1371 series Base Pole: 1" Diameter, 16 GA; Extension Pole: 1" Diameter, 18GA
Welding	All welding by heliarc method and exposed welds are ground smooth
Finish	Exposed surfaces are free from pit marks, weld seams, scale and are polished with #4 polish to a satin finish
Casters	See page 9
Notes	<ul style="list-style-type: none">MR conditional IV stand option available where indicated

Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs at info@blickman.com.

Note for: *IV Stands* Parts shipped separately to avoid damage and lower costs. Shipping weight with an asterisk is the combined shipping weight of separate parts. Some assembly required.

THE BLICKMAN DIFFERENCE

IV Irrigation Safety Systems

- Our IV Irrigation Safety Systems (IVTs) offer superior safety, efficiency, and financial value for fluid intensive IV irrigation procedures, such as arthroscopic, hysteroscopic, endoscopic, TURP, etc.
- They maintain continuous fluid delivery by offering infinitely adjustable heights for up to eight 3 liter or 5 liter IV bags on four independent channels.
- Their easy, one handed operation allows nurses to attach IV bags at comfortable heights, starting at 5 ft., and then raise bags up to 8 ft. (or to 9 ft. with custom extender) simply by pushing down on “squeeze/release” handles, an ergonomic solution that meets OSHA standards and greatly reduces the threat of back injuries.

Key Product Features

Up to eight 3 liter or 5 liter IV bags on four independent channels


High quality, easy maintenance, light weight aluminum

Locking casters for safe mobility and 360° spin in tight spaces

Push down “Squeeze/Release” handles to raise IV bags

Indefinite adjustable heights from 5 ft. to 8 ft. (9 ft. with custom extender)

26 1/2” diameter, inverted, 5 legged base resistant to tipping

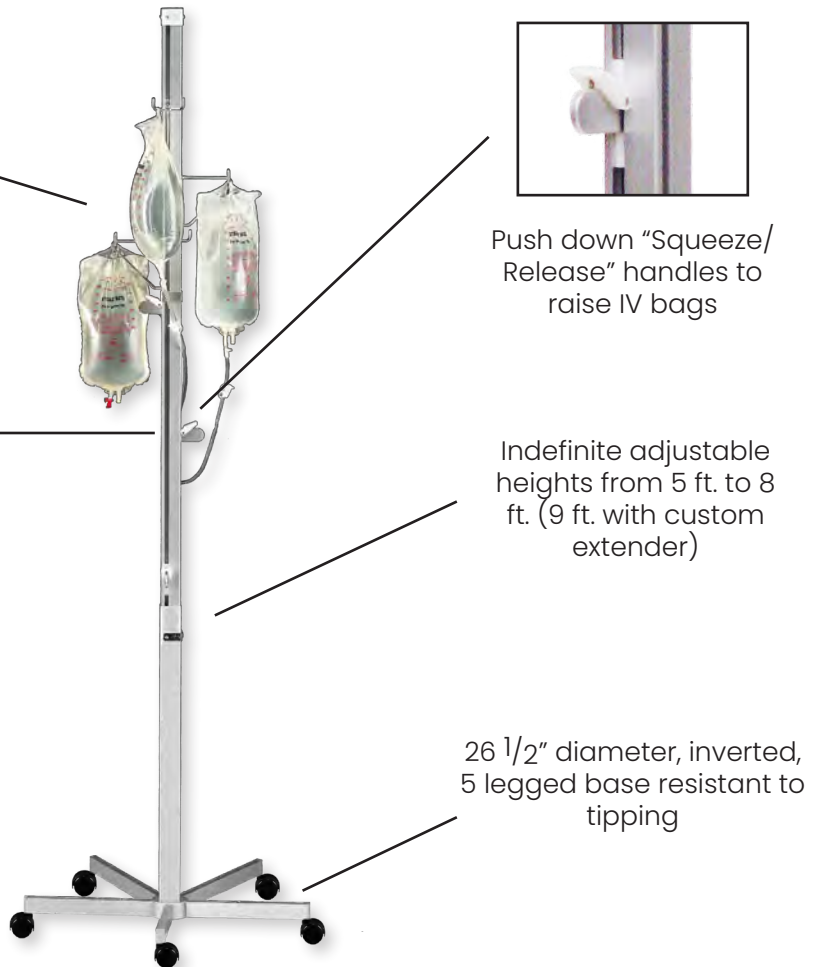

WHY CHOOSE OUR IVT?



Safe:

- This ergonomic solution meets OSHA standards and greatly reduces the threat of back injuries with its one-handed operation which allows nurses to attach IV bags at comfortable heights, simply by pushing down on “squeeze/release” handles, an
- Hard anodized aluminum construction with stain resistant finishes makes it so that no surfaces can chip, flake, or rust to introduce particle matter into a sterile or surgical environment.
- The lowest positions are cushioned by shock absorbers, thus sparing the common banging known to occur with the standard IV poles, which results in tearing of bag tabs and spilling remaining fluids onto the O.R. floor.
- The steady, five-legged bases resist potential tipping hazards, eliminating the need for “sand bagging”.

Efficient:

- There are no electrical aspects that require overnight charging, so they are always ready, whenever you are ready.
- The ability to stagger heights of up to eight 3 liter or 5 liter bags ensures continuous flow of liquids, as it allows you to change an IV bag, while the other bags are still irrigating.
- Easy to clean and maintain
- Compact footprint frees up valuable space on the O.R. floor.

Affordable:

- Our IV Irrigation systems are value leaders among a few comparable models, saving you up to **\$50/procedure** on IV tubes vs. those required with far more expensive, pump driven systems.

Irrigation Systems


IVTs

Description: We offer two types of IVTs, one with 2 hooks, which holds up to 4 IV bags and one with 4 hooks, which holds up to 8 IV bags. Hooks can raise and lower from 67" to 81". They travel on black plastic twin wheel hooded swivel casters. Both IVTs come with 26 1/2" diameter inverted, 5 leg base.

Our IVT stands have the same dimensions:

Width: 12 3/4" x Height: 89" x Depth: 12 3/4"

Note: these dimensions do not include the base

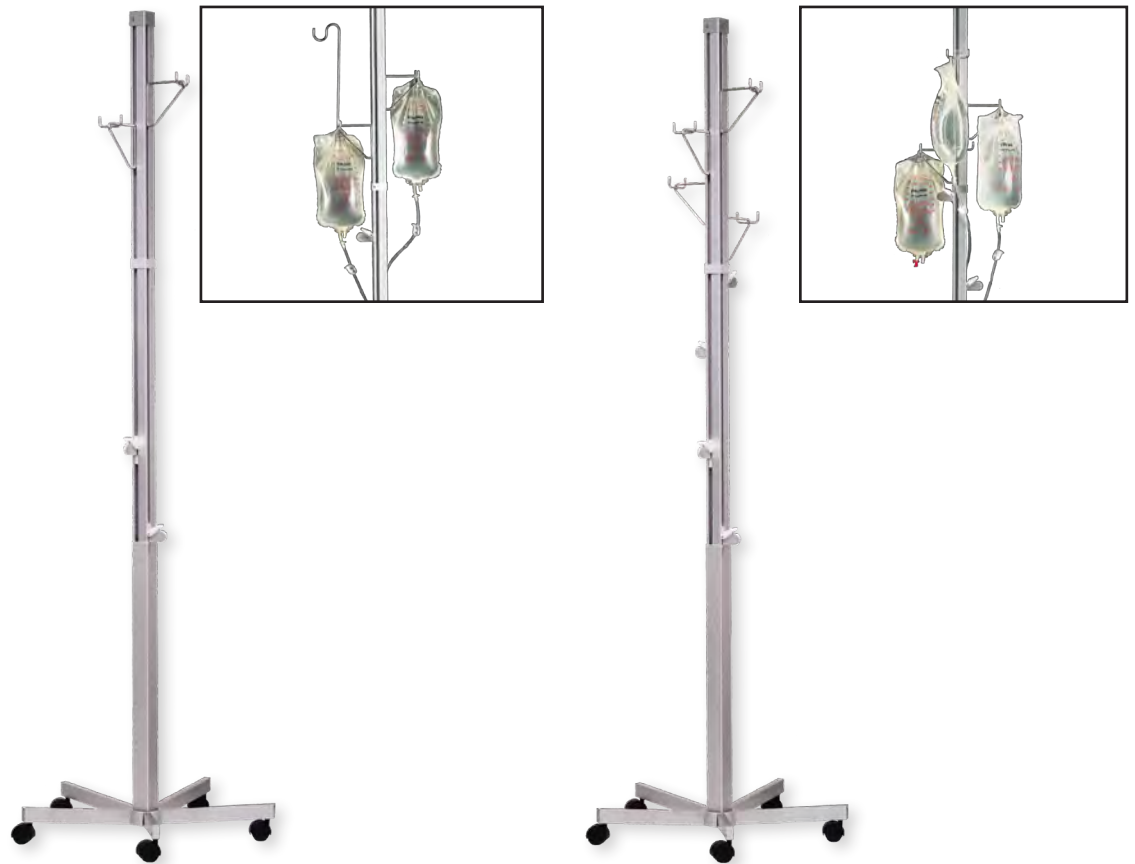

Name	IVT 2 Channels	IVT 4 Channels
Item #	IVT-2	IVT-4
Model #	0541110200	0541110100
Shipping Weight	85 lbs.*	85 lbs.*

For Custom Extender

Please email us at info@blickman.com

Irrigation Systems

Construction Specifications

Construction Specifications

Universal Dimensions Tower: Width: 12 3/4" x Height: 89" x Depth: 12 3/4"
Base Diameter: 26 1/2"
Height Range of Hooks: 67" - 81"

Materials Hard anodized aluminum

Gauge N/A

Welding N/A

Finish Stain resistant finishes

Casters 2" black plastic, twin wheel, hooded swivel casters

Note for: Parts shipped separately to avoid damage and lower costs.
IV Stands Some assembly required.

For Custom Extender

Please email us at info@blickman.com

Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs at info@blickman.com.

Toll Free: (800) 247-5070

Local: (973) 330-0557

BLICKMAN

www.blickman.com
info@blickman.com

1099 Wall Street West, Suite 250A, Lyndhurst, NJ 07071